

# Micropipet Protocol Card - **FINNPIPETTE**

**NEVER** crank the micropipets above or below their ranges!!

Range:  
**1.0-10  $\mu$ l**



Example  
volumes:




**1.2  $\mu$ l**



**6.0  $\mu$ l**

Range:  
**2-20  $\mu$ l**

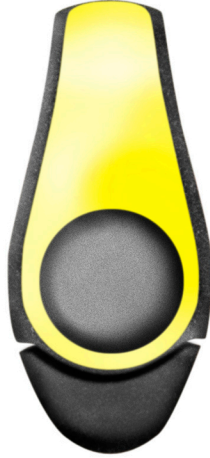

Example  
volumes:



**12.5  $\mu$ l**



**9.75  $\mu$ l**

Range:  
**20-200  $\mu$ l**



Example  
volumes:



**25.0  $\mu$ l**



**183.0  $\mu$ l**

Range:  
**200-1000  $\mu$ l**



Example  
volumes:




**1000.0  $\mu$ l**




**468.0  $\mu$ l**

Always hold the micropipet in a **vertical position** when there is fluid in the tip.

## **Tip Selection Guide**

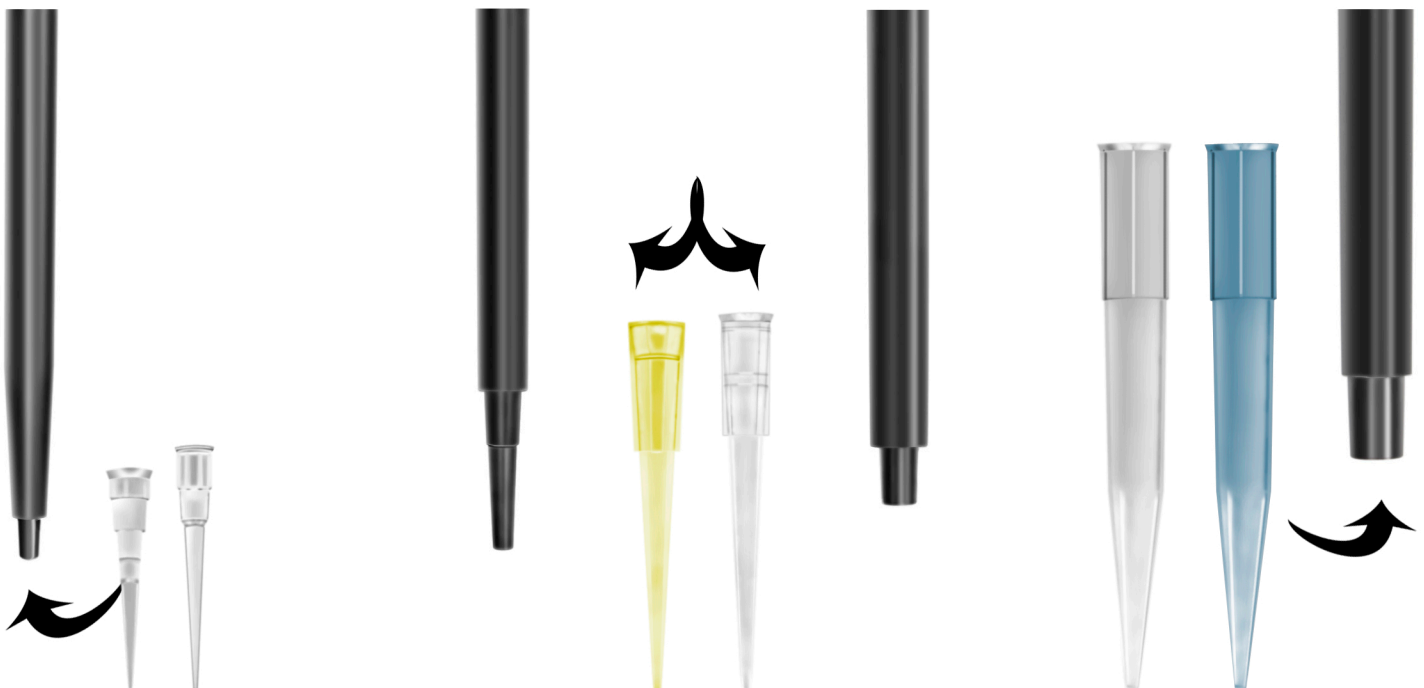

# Micropipet Protocol Card - **FINNPIPETTE**

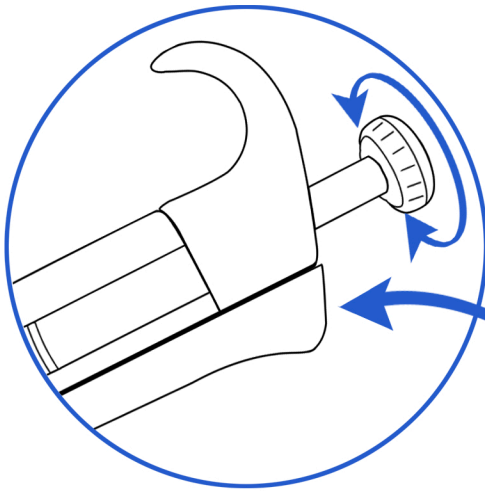

Set the volume by turning the Colored Knob

Tip Ejector Knob

Push down to first stop

Place tip into sample. Release thumb to draw up

Move to new microtube

Push down to first stop

Continue down to second stop

Eject Tip

